



# VERO CARE CLINER COLLECTION<sup>TM</sup>

OWNER'S OPERATING &  
MAINTENANCE MANUAL



Model 6Y48-39-00-00-00-D6  
Shown

PLEASE DO NOT DISCARD THIS MANUAL! KEEP FOR FUTURE REFERENCE AND TRAINING

# TABLE OF CONTENTS

<b>TABLE OF CONTENTS.....</b>	<b>2</b>
<b>INTENDED USE STATEMENT.....</b>	<b>3</b>
<b>TRANSPORTATION, STORAGE, HANDLING AND DISPOSAL INSTRUCTIONS.....</b>	<b>3</b>
<b>IMPORTANT SAFETY INFORMATION.....</b>	<b>3-6</b>
<b>ASSEMBLY INSTRUCTIONS.....</b>	<b>7-8</b>
<b>ATTENDANT OPERATING INSTRUCTIONS.....</b>	<b>9-14</b>
<b>OCCUPANT OPERATING INSTRUCTIONS.....</b>	<b>15-16</b>
<b>GENERAL MAINTENANCE.....</b>	<b>17</b>
<b>SPECIFICATIONS.....</b>	<b>18</b>
<b>VERO CARE CLINER OPTIONS.....</b>	<b>19-20</b>
<b>CONTACT &amp; WARRANTY INFORMATION.....</b>	<b>BACK COVER</b>

# INTENDED USE STATEMENT

**Intended Use:** Vero is a recliner for use in medical settings. It is designed to provide comfort for patients seated up to 8 hours at a time while receiving treatments such as Dialysis, Infusion or recovery from surgical procedures where positional adjustments are beneficial. The recliner is also suitable for general use in many medical settings such as, but not limited to hospital room chair, waiting room chair, convalescent home chair, etc.

## TRANSPORTATION, STORAGE, HANDLING AND DISPOSAL INSTRUCTIONS

The product should be **transported** in factory packaging, inside an appropriate medium for the destination, i.e., air/sea cargo containers. The product should be **stored** in an environment that will inhibit rust/mold formation. Winco also recommends not leaving the product in the factory packing in excess of 3 months. The product should always be **handled** in a manner consistent with the user instructions, and in a manner to prevent contamination after each use. The product has many recyclable components and to the extent practical, all effort should be used to recycle responsibly. Otherwise the product components shall be **disposed** of in accordance with local statutes.

***Winco assumes no responsibility for damage or injury caused by the improper assembly, installation, use or maintenance of this product.***

### **IMPORTANT!**


PLEASE READ THIS ENTIRE MANUAL BEFORE USING RECLINER. DO NOT INSTALL, MAINTAIN OR OPERATE THIS EQUIPMENT WITHOUT READING AND FOLLOWING THIS MANUAL. OTHERWISE, SERIOUS INJURY AND/OR DAMAGE TO PRODUCT MAY OCCUR.

THE INFORMATION CONTAINED IN THIS MANUAL IS SUBJECT TO CHANGE WITHOUT NOTICE. PHOTOS ARE REPRESENTATIVE OF THE PRODUCT AND MAY VARY SLIGHTLY FROM ACTUAL PRODUCTION MODELS.

### **SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE!**

NO PART OF THIS MANUAL MAY BE DUPLICATED IN ANY FORM WITHOUT THE PRIOR WRITTEN CONSENT OF WINCO MFG, LLC. UNAUTHORIZED DUPLICATION MAY RESULT IN CIVIL PROSECUTION TO THE MAXIMUM EXTENT ALLOWED BY LAW.

***Products with the following markings are intended for commercial use only.***

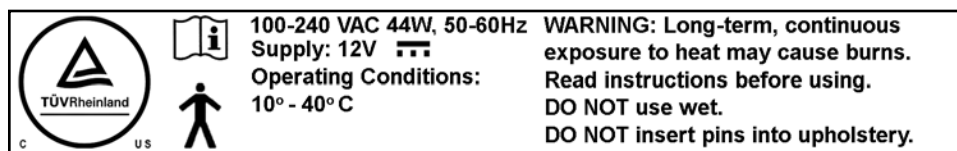

### **SYMBOLS**



FOLLOW INSTRUCTIONS



DC VOLTAGE



TYPE B APPLIED PART



GENERAL WARNING SIGN



PINCH POINT  
WARNING SIGN



TRENDELENBURG



SYMBOL FOR MANUFACTURER

# IMPORTANT - PLEASE READ



**FOR PRODUCTS WITH POWERED OPTIONS (Heat, Massage, etc...)**

**IMPORTANT SAFETY INSTRUCTIONS - SAVE THESE INSTRUCTIONS**

***(Read all instructions before using this product.)***

***When using an electrical appliance, basic precautions should always be followed, including the following:***

**DANGER** - To reduce the risk of electric shock always unplug this product from the electrical outlet before cleaning, maintenance or servicing.

**WARNING:**

1. Unplug this product from the electrical outlet when not in use.
2. Supervision should be provided when using this product.
3. Use this product only for its intended use as described in this manual. See page 3 for intended use.
4. Do not use attachments not recommended by the manufacturer.
5. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water, contact an authorized Winco Representative for instructions.
6. DO NOT carry or pull this product by power supply cord or use the cord as a handle.
7. Keep the cord away from heated surfaces.
8. Do not use outdoors.
9. To disconnect, turn all controls to the "off" position, then remove plug from outlet. Pull plug only to remove cord.
10. Clean switch with a damp cloth only – Warranty will be VOID if liquid is introduced into switch mechanism.
11. Extension cords are not recommended.
12. DO NOT roll recliner or other equipment over cord.
13. This product is for use on a nominal 120-volt circuit and has a three-pin plug. If the plug does not fit into the outlet, contact a certified electrician to install the proper outlet. DO NOT MODIFY THE PLUG OR USE AN ADAPTER WITH THIS PRODUCT.
14. For an added level of safety; Winco recommends the use of a 120 volt GFCI (ground fault circuit interrupter) outlet.

**DANGER:**

- i. Never insert pins, or attach other metallic fasteners into/onto any part of this chair.
- ii. Do not use this product if the covering shows signs of deterioration, such as checking, blistering, or cracking.
- iii. KEEP DRY - Do not operate in a wet or moist condition or environment.


**WARNING:**

- i. Long term, continuous exposure to heat may cause burns. It is recommended that skin in contact with heated area be checked for redness and blistering during long term usage. Do not use on an infant, invalid or incapacitated person. Do not use on insensitive skin or on a person with poor blood circulation. If you are uncertain if you should use this product, please consult your physician about the use of heat and massage.
- ii. Do Not Crush/Pinch heating elements or wiring.

# IMPORTANT - PLEASE READ


For Heat & Heat and Massage Chairs

This equipment is marked with the industry recognized symbol

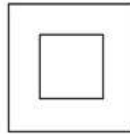

indicating it is a "Class II Double-Insulated" device as defined by IEC 60601-1. As such, it does NOT require a three conductor power cord (i.e. one with a ground connection) for safety purposes. While the equipment can be safely operated with no ground connection, a three conductor power cord is provided to satisfy customer expectations for typical equipment configuration. Devices qualifying as double insulated require different testing & examination procedures than "Class I Grounded" devices. If your facility requires the product to be tested, be sure to follow testing & examination procedures for "Class II Double-Insulated" devices.

ISOLATION DIAGRAM FOR HEAT & MASSAGE SYSTEM

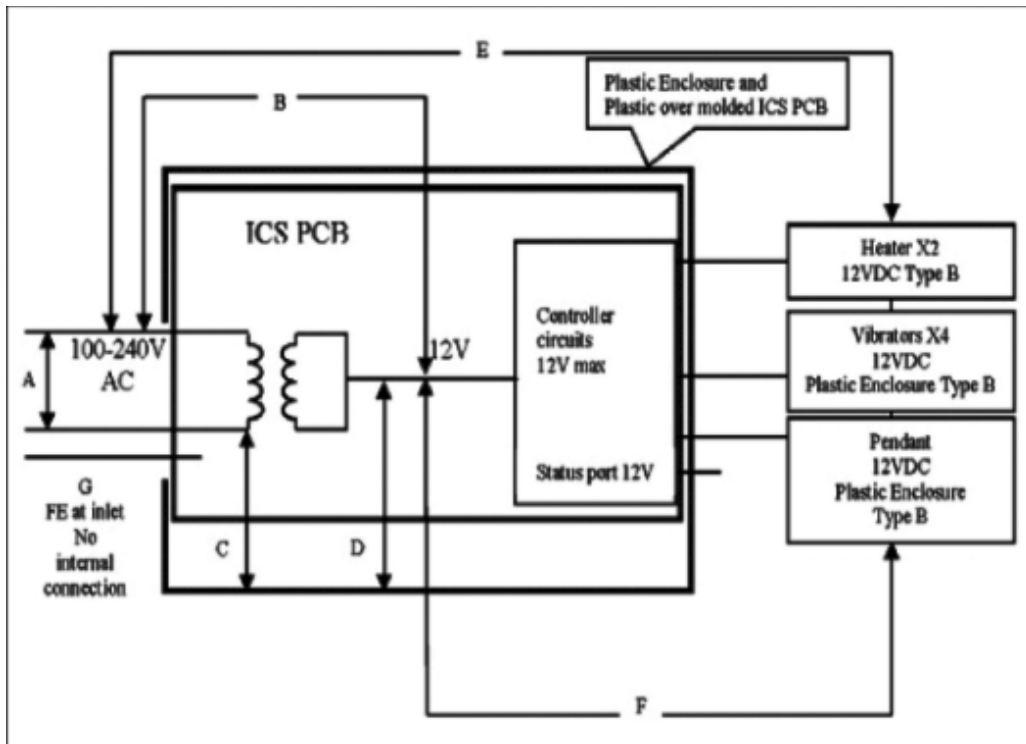

Area Insulation Type

- A Basic Operational
- B Double Insulated/Reinforced
- C Double Insulated/Reinforced
- D Double Insulated/Reinforced
- E Double Insulated/Reinforced
- F Double Insulated/Reinforced
- G Inlet FE pin only connects to inlet & pad


# IMPORTANT - PLEASE READ



## **BASIC SAFETY PRECAUTIONS:**

1. READ AND FOLLOW ALL DIRECTIONS.
2. CAUTION: Use slide-out foot support with staff assistance only.
3. DO NOT put hands, feet, or clothing into any openings when changing positions on recliner. Attendant MUST confirm that the occupants arms, legs, hands and feet are clear while changing recliner positions or SERIOUS INJURY MAY OCCUR.
4. Remove or move chair accessories out of the way of patients legs (ie. over-lap table & pivot table) when reclining/ changing position of patient.
5. LOCK casters at all times when chair is stationary. UNLOCK casters when moving/transporting chair.
6. STAY CLEAR of the recliner mechanisms when operating chair. Keep children away from extended foot support (or other similiar parts).
7. DO NOT use the recliner if the Trendelenburg Damper/Gas Spring is removed or damaged.
8. DO NOT use recliner for transporting in or with ANY type of vehicle or trailer. Winco recliners have not been tested or approved for use by an occupant in any type of vehicle or trailer.
9. NEVER use the chair arms or backrest as a seat; SERIOUS INJURY OR DAMAGE MAY OCCUR.
10. Periodically check the tightness of all nuts, bolts and screws.
11. Immediately REMOVE FROM SERVICE any recliner with broken recline mechanisms, missing parts, torn upholstery, and/or other mechanical or visible damage.
12. Transport user with chair arms in the open position only in emergency situations when needed to get chair through narrow exit doors.
13. USE ONLY WINCO AUTHORIZED REPLACEMENT PARTS.
14. NEVER EXCEED the recommended weight capacity.

**WARNING:** USE CAUTION WHEN CLEANING NEAR MECHANISMS. KEEP FINGERS AWAY FROM PINCH POINTS.

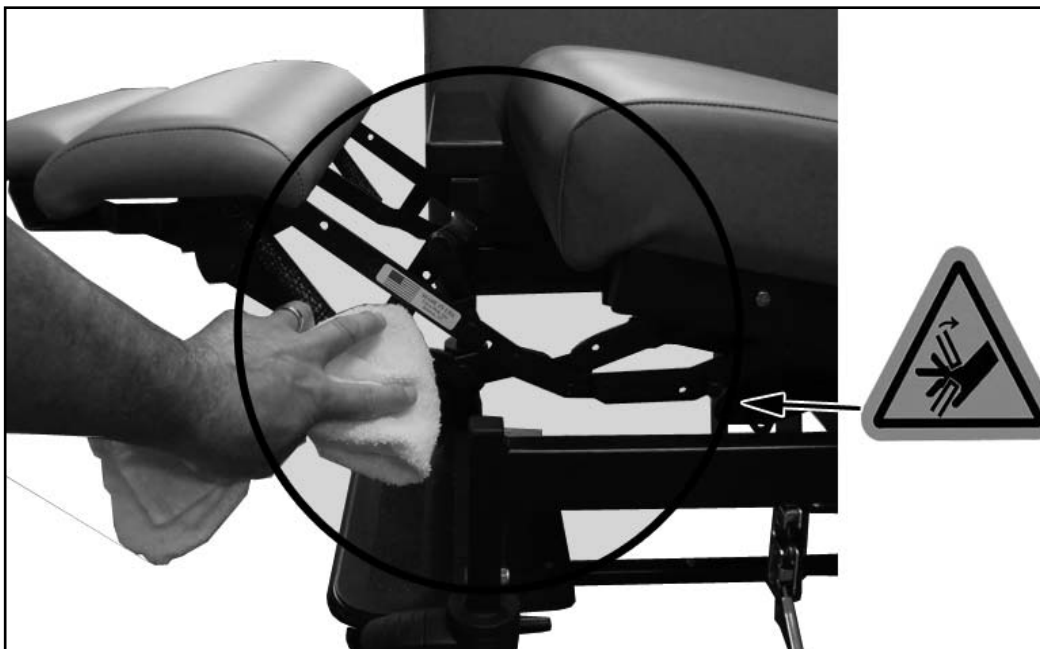

# ASSEMBLY INSTRUCTIONS

## PREPARATION:

1. Carefully examine your product for any damage. Inspect all components.  
IF DAMAGE IS EVIDENT, REFUSE SHIPMENT, & CONTACT FREIGHT CARRIER OR WINCO IMMEDIATELY.
2. Remove all packaging material that was used for shipping.
3. Carefully remove the carton from the product.
4. Save all boxes and packaging material until after you have assembled your product and have verified that all components are functioning properly. These material are required if it becomes necessary to return the product.

## BACK ASSEMBLY INSTRUCTIONS:

1. Lock ALL (4) casters on chair base. (See Lock/Unlock Casters pg.8)

### IMPORTANT:

- Keep objects and persons clear of the chair before performing step 2.

2. Recline mechanism should be in the upright position. Remove any packaging that may be on the mechanism that might prevent installation of back.



3. Hold the rear cover, on the lower rear of the back, up exposing the mounting brackets (A). Lift the upholstered back by its sides while aligning the mounting brackets with the recline mechanism (B).

### WARNING:

- Keep hands clear of connection points or injury may occur

### IMPORTANT:


- Winco recommends two people be used for step 4
- Be careful not to tear or puncture vinyl.

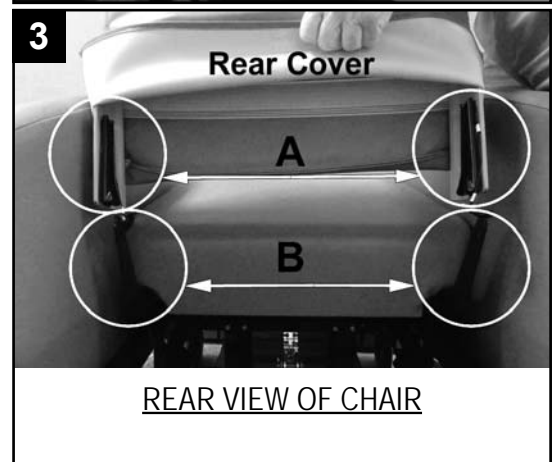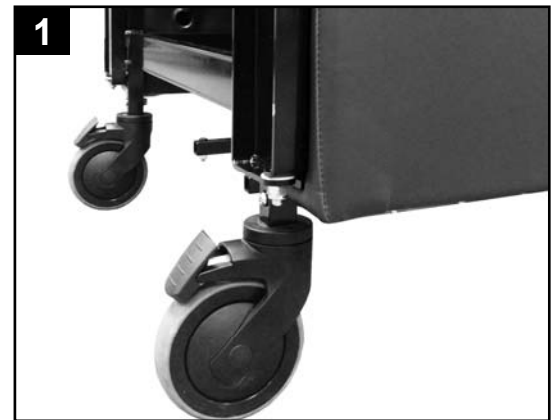


# ASSEMBLY INSTRUCTIONS

## BACK ASSEMBLY INSTRUCTIONS:

5. Slide the back onto the recline mechanism until the middle hole of the mounting bracket (A) is aligned with the middle hole of the recline mechanism (B).

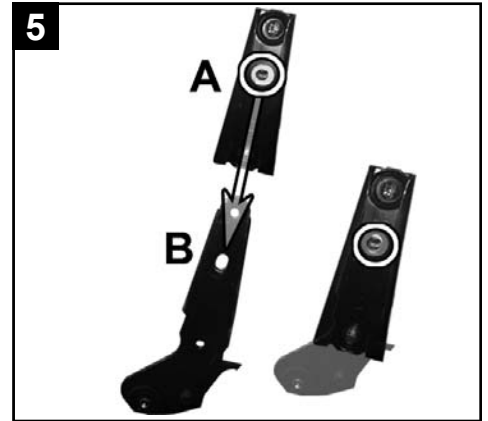

6. Install knobs (tethered to back) and fully tighten to a secure fit.

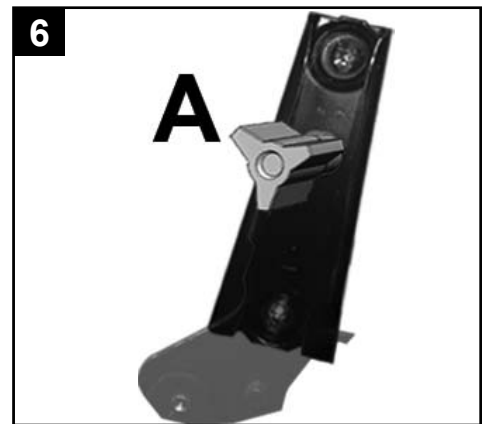

### IMPORTANT:

- Test to ensure the back has been properly installed by pulling up on the back. The back should remain securely in place.



# ATTENDANT OPERATING INSTRUCTIONS


## CASTER OPERATION: Lock/Unlock

This chair offers multiple types of casters: rigid and total-lock (5" black nylon), wheel-lock (3" black nylon) and pedestal feet.

### Rigid and Total-Lock:

- The **front** two casters are rigid casters. They are NOT designed to swivel in order to facilitate stable straight-line steering. The front caster locking tabs are always accessible. When locking tab is engaged, caster wheel will NOT roll. The **rear** two casters are total-lock casters. They are designed to swivel 360°. When locking tab is engaged, caster wheels will NOT roll and caster will NOT swivel.

1. TO LOCK CASTER: Press down on the wheel tab (A). Tab will lock into position shown.
2. TO UNLOCK CASTER: Press down on the TOP of the wheel tab (B). Tab will return to position shown.

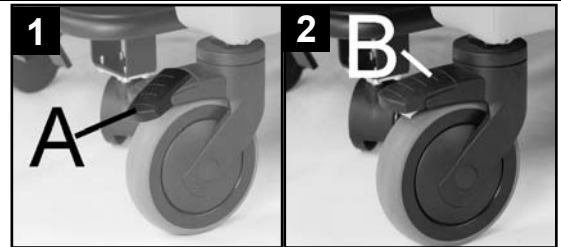

### Total-Lock:

- Total-Lock casters lock the wheel and swivel mechanism when pressing down on the tab.




**WARNING: LOCK CASTERS AT ALL TIMES,  
EXCEPT WHEN MOVING/TRANSPORTING CHAIR**

### Wheel-Lock casters (3" Black casters only):

- Wheel-Lock casters (3" Black casters only) will lock wheel but does NOT prevent wheel from swiveling.
1. TO LOCK CASTER: Press down on the wheel tab (A). Tab will lock into position shown.
  2. TO UNLOCK CASTER: Press down on the TOP of the wheel tab (B). Tab will return to position shown.

**Note:** All 3" Black casters are wheel locking ONLY.  
TOTAL and DIRECTIONAL lock features are not available.

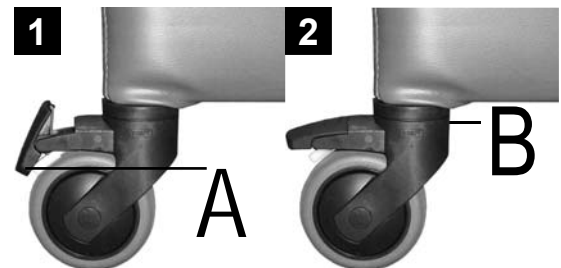

**WARNING: ALWAYS LOCK THE CASTERS WHEN THE OCCUPANT ENTERS OR EXITS THE CHAIR! UNLOCK CASTERS WHEN TRANSPORTING CHAIR!!**

### Pedestal Feet:

- Pedestal feet are designed for chairs that are intended to remain in a stationary position.

**Note:** The Foot Plate option is not available with Pedestal Feet or 3" Black Casters.

# ATTENDANT OPERATING INSTRUCTIONS

## RECLINE CHAIR FROM ATTENDANT POSITION:

Attendant should be positioned on the RIGHT or LEFT side of chair

### CHAIRS WITH MANUAL ADJUSTMENT

1. While occupant is seated in the chair, place one hand on the arm for support and push back on the back to recline.

### CHAIRS WITH INFINITE GAS SPRING ADJUSTMENT OPTION


(These models have the recline lever (A))

1. To operate the infinite back positioning use one hand to operate the recline lever (A) and the other hand to push on the upholstered back (B) until desired position or maximum recline has been reached.

Release recline lever (A) to maintain desired position.

#### IMPORTANT:

- Do NOT force the back of the chair up or down without operating the recline lever, or damage to the chair may occur.

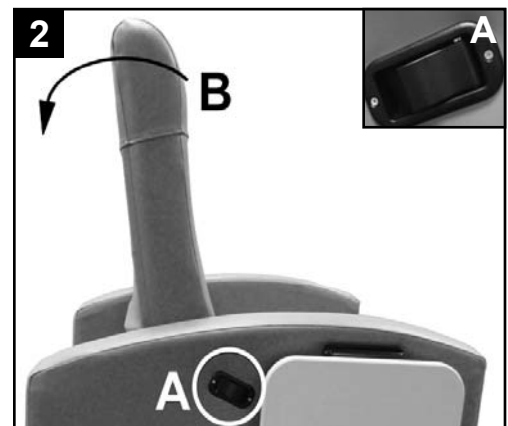

## RETURN CHAIR UPRIGHT FROM ATTENDANT POSITION:

Attendant should be positioned on the RIGHT or LEFT side of chair

### CHAIRS WITH MANUAL ADJUSTMENT

1. Pull up on the back to return to upright position.
2. Attendant can use one hand to push down on the legrest (C) until it returns to its stored position. (**Warning - Keep hands clear of mechanism to avoid injury**)

### CHAIRS WITH INFINITE GAS SPRING ADJUSTMENT OPTION

(These models have the recline lever (A))


1. While occupant is reclined in the chair, attendant should use one hand to operate the recline lever (A), while using the other hand to pull up on back (B).
2. Attendant can use one hand to push down on the legrest (C) until it returns to its stored position. (**Warning - Keep hands clear of mechanism to avoid injury**)

Release the recline lever (A) once the back has fully returned to the upright position.

**NOTE:** The recline lever is NOT required to return legrest to its stored position.

#### WARNING:

- Attendant should NOT use their foot to return the legrest to its stored position. Damage to mechanism may occur.

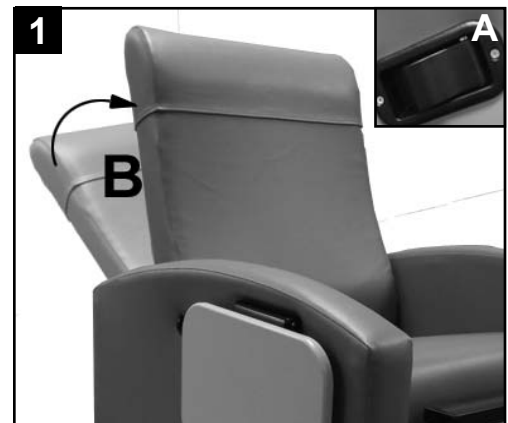

# ATTENDANT OPERATING INSTRUCTIONS

## PLACING CHAIR INTO TRENDELENBURG: (*Attendant ONLY*) (*OPTIONAL FEATURE*)

Attendant should be positioned on the RIGHT or LEFT side of chair

1. While occupant is in the chair, place the chair into the fully reclined position. (SEE "Recline Chair From Attendant Position" pg.9)



2. Once chair is fully reclined, attendant can push down on the rear of Trendelenburg lever (A) with their foot, allowing it to move down.
3. Hold Trendelenburg lever in position & push down on the chair's back (B). This will bring the occupant into the Trendelenburg position.

Release the Trendelenburg lever.

**NOTE:** The chair does not lock in the Trendelenburg position.

**IMPORTANT:**

- Chair must be in the fully reclined position BEFORE occupant can be placed into Trendelenburg.




## FROM TRENDELENBURG TO RECLINE: (*Attendant ONLY*)

Attendant should be positioned on the RIGHT or LEFT side of chair

1. To return chair to the recline position, push down on seat (recommended) or lift up on chair back.

**NOTE:** Pushing down on the seat or legrest (B) is recommended for maximum leverage.

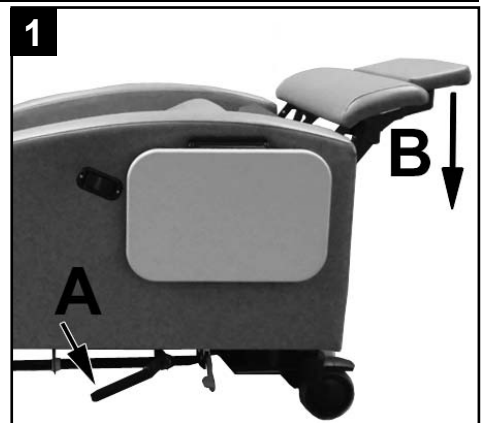

# ATTENDANT OPERATING INSTRUCTIONS


## SWING-ARM OPERATION: (*Attendant ONLY*) (*OPTIONAL FEATURE*)

Attendant should be positioned on the RIGHT or LEFT side of chair

1. TO OPEN: Press down on red latch handle to unlock Swing-Arm.

**WARNING:**

- Do NOT leave occupant unattended with either or both swing-arms unlatched.
- Do NOT transport chair with swing-arms unlatched unless it is an Emergency situation.

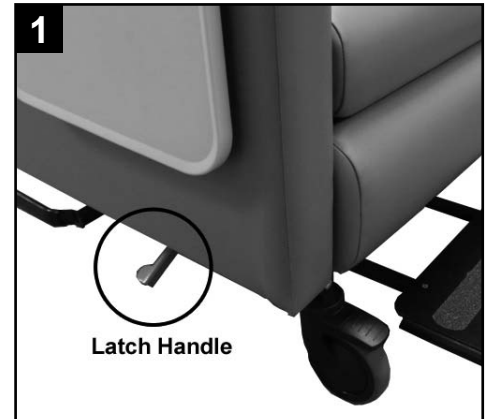

2. TO CLOSE: Push Swing-Arm firmly against frame.  
Latch will automatically engage with an audible "click".

**WARNING:**

- Make sure swing-arm latches are properly engaged when swing-arms are in the closed (latched) position.

Test by pulling on arms. The arms should not swing open if properly engaged.



### 180° SWING-ARM OPTION:

This chair option provides swing-arms that pivot up to 180°. When both arms are pivoted fully back, chair width is reduced to allow the chair to be more easily transported through narrow doorways. In addition, swing-arm can be pivoted back in order to get a wheel chair or gurney alongside the chair. (See pg. 12)

3. Unlock both swing-arms and rotate them towards the back of the chair.
4. Stand behind the chair and use the back or push handle (if equipped) to maneuver the chair through the doorway.

**IMPORTANT:**

- Transport occupant with chair arms in the open position in emergency situations only. This should only be performed by an attendant.

**WARNING:**

- Do NOT recline chair with either/both swing-arms in the 180° position.
- Ensure that both swing-arms are in the closed (latched) position at all times, unless the patient is being transferred or directly supervised.



# ATTENDANT OPERATING INSTRUCTIONS

## TRANSFER OF OCCUPANT TO CHAIR (Swing-Arm Option Only):

Attendant should be positioned on the RIGHT or LEFT side of chair

**WARNING: ALWAYS LOCK THE CASTERS WHEN THE OCCUPANT ENTERS OR EXITS THE CHAIR!**

1. Lock all four casters.
2. Place chair into seated or reclined position relative to the position from which occupant is being transferred. (See pg.9)
3. Unlatch one swing-arm and place it into the 180° position.
4. Position gurney or wheel chair next to the chair on the side with the open swing-arm.
5. Transfer patient to chair following facility safe patient handling policies.
6. Close the swing-arm and ensure latches engage properly.  
(See pg. 11-Item 2.)

**WARNING:**

- Do NOT leave swing-arm open after occupant is transferred to chair.
- When closing swing-arm, make sure occupants extremities are clear of swing-arm.
- When closing swing-arm, make sure any foreign objects are clear of swing-arm travel and latch.


## FOOTPLATE OPERATION: *(OPTIONAL FEATURE)*

**IMPORTANT:**

- Footplate MUST be used with staff assistance.
  - The footplate for the VERO model is rated for a MAXIMUM weight capacity of 350 lbs. (158.8 kg.)
  - The footplate for the VERO XL model is rated for a MAXIMUM weight capacity of 500 lbs (226.8 kg.)
1. TO EXTEND: Using your foot, pull out footplate as far as it will travel. Footplate can now be used by the occupant.
  2. TO STORE: Using your foot, push in on the footplate until it is returned to the fully stored position

**WARNING:**

- Footplate can be a tripping hazard if not returned to its stored position after use.
- ALWAYS store footplate when not in use.

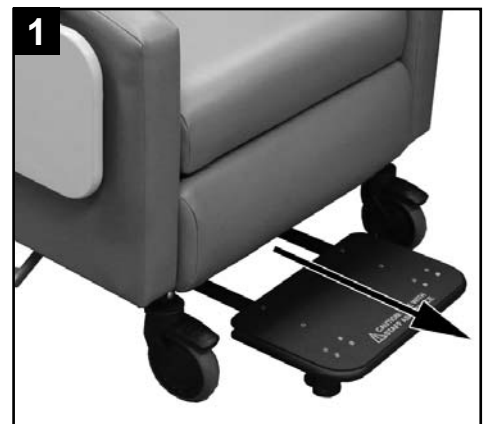


# ATTENDANT OPERATING INSTRUCTIONS

## SIDE AND PIVOT TABLE OPERATION: *(OPTIONAL FEATURE)*


### **PREMIUM STYLE SIDE TABLE:**

1. **TO RAISE TABLE:** Grasp center of table and gently lift into place. Side-table will lock into place with an audible "click" once fully extended.
2. **TO STORE:** Lift up on the lever located centrally under the side-table and lower to its stored position.

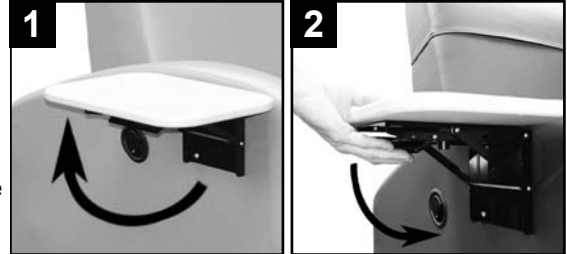

### **STANDARD SIDE TABLE:**


1. **TO RAISE TABLE:** Grasp center of table and gently raise into place. Table will automatically lock into place when released.
2. **TO LOWER TABLE:** Grasp table at the center, or on both sides, raise and pivot down in one movement.

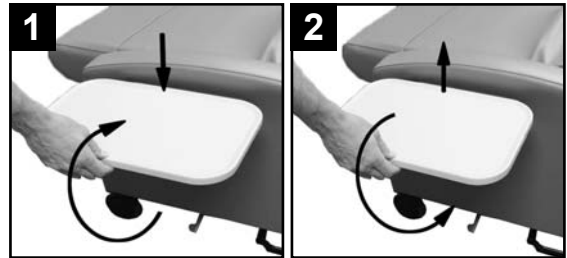

### **PIVOT TABLE:**

#### **TO POSITION THE PIVOT TABLE:**

1. Loosen the adjustment knob (A) and swing away or place in the "lap" position.
2. Tighten adjustment knob (A) to secure position.

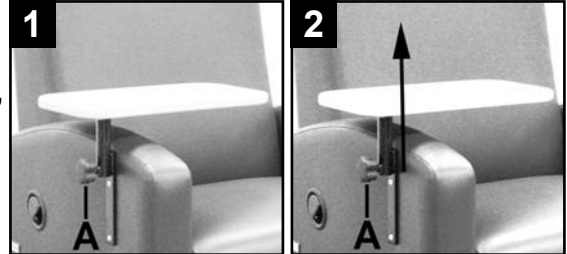

#### **TO REMOVE THE PIVOT TABLE:**

3. Loosen the adjustment knob (A) 2 full turns and lift up on pivot-table.

#### WARNING:

- DO NOT lean or sit on any Vero tables.
- DO NOT use table for support or assistance when exiting or entering the chair.
- Keep hands and objects clear of the moving parts of the table mechanism.

#### IMPORTANT:

- Side and Pivot tables have been rated for a MAXIMUM weight capacity of 20 lbs. (9.0 kg)

# OCCUPANT OPERATING INSTRUCTIONS

## RECLINE WHILE IN THE SEATED POSITION:

### CHAIRS WITH MANUAL ADJUSTMENT

1. While seated in the chair, place hands on each chair arm and push back while applying pressure to the back with shoulders to activate the mechanism into the desired position.



### CHAIRS WITH INFINITE GAS SPRING ADJUSTMENT OPTION


(These models have the recline lever (A))

1. To operate the infinite back positioning, pull back on the recline lever and adjust the upholstered back angle by pushing back with your body until desired angle has been reached.

Release recline lever (A) when you have achieved the desired angle to lock the back into position.

#### WARNING:

- DO NOT stand or sit on the legrest
- DO NOT place hands, legs, or feet under seat, mid-ottoman, or legrest.
- Stay clear of recline mechanism when operating chair.

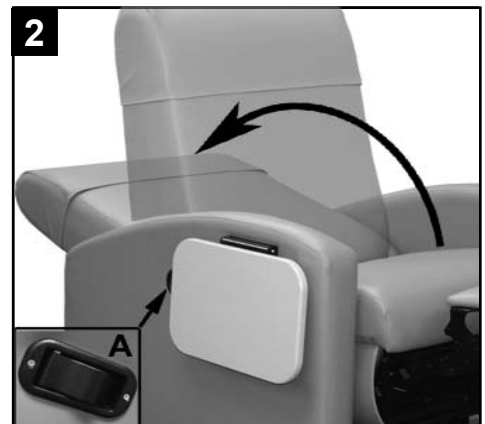


## RETURN TO THE SEATED POSITION:

### CHAIRS WITH MANUAL ADJUSTMENT

1. To return the back to the full upright position, grasp chair arms with both hands and pull body forward. (Go to step 2)
2. Press down on the legrest (A) with legs/feet until it is stored.  
**WARNING: Make sure all objects, hands and feet are clear mechanisms before lowering legrest!**

**NOTE:** Use of the recline lever is NOT required to return the legrest to its stored position.

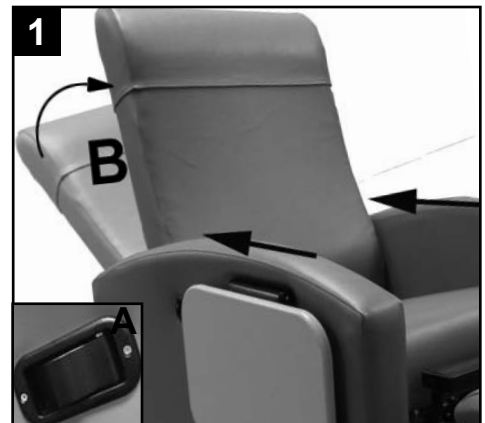

### CHAIRS WITH INFINITE GAS SPRING ADJUSTMENT OPTION

(These models have the recline lever (A))

1. To return the back to the full upright position, pull back on either of the recline lever (A) and lean body forward.

Release the recline lever once the back has fully returned to the upright position.

2. Press down on the legrest (A) with legs/feet until it is stored.  
**NOTE:** Use of the recline lever is NOT required to return the legrest to its stored position.

**(WARNING: Make sure all objects, hands and feet are clear mechanisms before lowering legrest!)**



# OCCUPANT OPERATING INSTRUCTIONS

## OPTIONAL HEAT & MASSAGE OPERATION:

Control is located on the inside of right arm of chair. (Right is determined - as if sitting in the chair)



**Heating:** (Heating elements are located in the seat & back of chair)

Push HEAT button - See indicator light on switch

**3 Lights** = High Heat / Bright Light

**2 Lights** = Medium Heat / Med Light

**1 Light** = Low Heat / Dim Light

**No Light** = HEAT OFF

Allow 5-8 minutes for back and seat to warm to desired setting.

Heat will automatically shut off after approximately 60 minutes. (Timer resets each time switch is pressed)

**Massage option:** (The massage option uses 4 separate motors located in the backrest of the chair. Two in the lower backrest, and two in the middle (upper) backrest.

Push MASSAGE button

**1X** = Lower Back / Low Intensity

**2X** = Lower Back / High Intensity

**3X** = Upper & Lower Back / Low Intensity

**4X** = Upper & Lower Back / High Intensity

**5X** = Alternating Upper & Lower Back / High Intensity

**6X** = MASSAGE OFF

Massage will automatically shut off after approximately 15 minutes. (Timer resets each time switch is pressed)

### CLEANING:

- With system disconnected from power source, use a DAMP CLOTH ONLY to clean switch, wiring or power pack (located under seat frame).

### TROUBLESHOOTING:

- If there is no heat or massage or flashing light on the switch - check power cord connections to wall outlet AND connection of power cord to power pack unit under chair.
- If switch displays a flashing light; Disconnect from wall, check all accessible wire connectors under chair for possible loose connections, then re-plug into wall and test again.
- If system still does not work contact Winco customer service: 1-800-237-3377.

### CORD STORAGE:

(Models with Heat-HT or Heat & Massage-HM)

To prevent tripping on cord, damage to cord, damage to recliner, power cord should be coiled and stored in the "Hook & Loop" fastener provided BEFORE transporting chair or when heat & massage is not in use.





# GENERAL MAINTENANCE

## **GAS CYLINDER OPERATION:**

### *DO NOT ATTEMPT TO ADJUST CYLINDER!*

The Gas Cylinder is preset during manufacturing and should not be adjusted without contacting Winco Customer Service.

- Adjusting the cylinder will not make the mechanism operate easier or more smoothly.
- If cylinder is removed and allowed to extend it will not be able to be reinstalled without special tools.

## **GENERAL CARE & CLEANING:**

- OCCASIONALLY lubricate pivot points on the recline mechanism with a lubricant approved by your facility.  
(If equipped) Periodically position the chair back into full recline position to keep the gas cylinder lubricated.
- Periodically inspect tightness of all screws, bolts, and nuts.

Winco furnishings are built to provide durable reliable service when properly cared for.

In general all of our products should be 1) Cleaned 2) Rinsed 3) Allowed to air dry.

It is important to note that disinfecting a product is not the same as cleaning the product. Disinfectants alone will not provide adequate cleaning, they do not have grease or oil cutting properties to remove grime or hair and skin oils. Over time, the appearance and feel and performance of your upholstery may diminish if not cleaned properly

### **CLEANING:**

Remove grime, hair and body oils; Your furnishings should be cleaned with mild soap & water solution and a damp cloth on a regular basis (especially where skin & hair make contact with upholstery). Avoid harsh detergents that could damage the upholstery or finish of your recliner. If the furnishing is disinfected with bleach, it MUST be wiped off using only clean water on a damp cloth and then be allowed to air dry. Failure to rinse upholstery with clean water can result in a build-up of residues over time that can dry out the upholstery and cause changes to the look and feel of the upholstery and eventually can lead to cracking. When cleaning or disinfecting the inside of a swing arm, it is important to allow the arm to air dry before closing the arm or wipe with clean dry cloth before closing. Otherwise, the moisture between the arm and the frame could cause the frame to rust.

# VERO CARE CLINER SPECIFICATIONS



Model 6Y48-39-00-00-00-D6  
Shown

VERO CARE CLINER - Specifications	VERO XL CARE CLINER - Specifications
Seat Height .....21" to Foot Plate.....18"	Seat Height .....21" to Foot Plate.....18"
Seat Width (arm to arm).....22.5"	Seat Width (arm to arm).....26.5"
Overall Height .....	Overall Height .....
5" Casters.....48.25"	5" Casters.....48.25"
3" Casters.....46"	3" Casters.....46"
Pedestal Feet.....44.75"	Pedestal Feet.....44.75"
Overall Width	Overall Width
(No Trays w/o Trend) .....29.5	(No Trays w/o Trend) .....34"
(No Trays w/Trend) .....30.25"	(No Trays w/Trend) .....34.25"
(D6 Trays Down) .....32"	(D6 Trays Down) .....36"
(D6 Trays Up) .....51.25"	(D6 Trays Up) .....55.25"
(DT Trays Down) .....35"	(DT Trays Down) .....39"
(DT Trays Up) .....55.5"	(DT Trays Up) .....59.5"
(Arms at 180 w/Trend) .....30.25"	(Arms at 180 w/Trend) .....34.25"
(Arms at 180 w/No Trend) .....28.5"	(Arms at 180 w/No Trend) .....32.5"
Back Height Above Seat.....29.25"	Back Height Above Seat.....29.25"
Seat Depth.....22"	Seat Depth.....22"
Weight .....131 lbs. <b>(Weight includes Footplate and Sidetables)</b>	Weight .....153 lbs. <b>(Weight includes Footplate and Sidetables)</b>
<b>Maximum Weight Capacity....350 lbs.</b>	<b>Maximum Weight Capacity....500 lbs.</b>

# VERO CARE CLINER OPTIONS

## Heat & Massage

HM .....Heat & Massage  
 HT .....Heat ONLY  
 FACTORY INSTALLED ONLY



## Premium Style Side-Table

Available on LEFT, RIGHT, or BOTH sides  
 (Also available with Wood Finish)



## Standard Side -Table

Available on LEFT, RIGHT, or BOTH sides  
 (Also available with Wood Finish)



## Pivot-Table

Available on LEFT, RIGHT, or BOTH sides  
 (NOT available with Wood Finish)



## Accessory Holder

Available on LEFT, RIGHT, or BOTH sides



## IV Pole & Attachment

IV..... IV Pole & Attachment  
 (Available on Left, Right or Both Sides)



## Headrest Cover

Option on Non-Trendelenburg Models



## Push Handle

Available for Vero and Vero XL Models



## Foot Plate

Available with 5" Casters Models Only



# VERO CARE CLINER OPTIONS

## Pedestal Feet

Not Compatible with the Foot Plate Option



## 3" Caster

Not Compatible with the Foot Plate Option



## Urethane Arm



## Arm Rest Covers


Left and right available for all table and option configurations.



# VERO CARE CLINER ARMREST COVER INSTRUCTION

## Installing Armrest Covers With Side Tables:

1. Raise the side table by grasping the center of the table and gently lift into place.

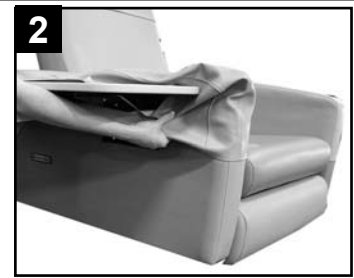


2. Locate large square opening in cover and position over the side table. Stretch the opening over the side table and slide over side table.

3. Fit the front portion of the Armrest cover over the front of the arm. Pull down tight for best fit.



4. Fit the rear portion of the Armrest cover over the rear of the arm.


5. For Premium Table, you will need to release the table slightly to allow space to tuck the apron into the gap. Standard Table will have enough room in the upright position.

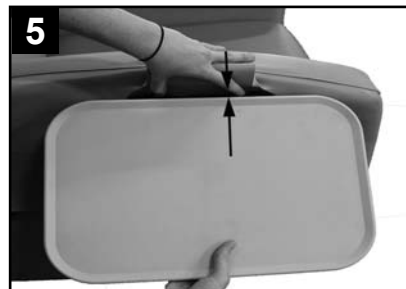

6. Pull the apron through the gap at the back of the side table completely. (Premium side table pictured) Smooth wrinkles.

## Installing Armrest Covers With Pivot Tables:


1. Remove the adjustment knob and the Pivot table completely. **You must remove the knob!**

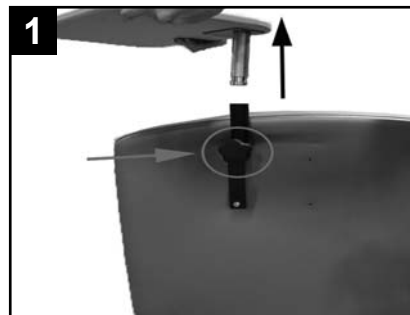

2. Pull the arm rest cover down over the chair arm. Ensure pivot table bracket is inserted through the hole in the arm rest cover.



3. Reinstall pivot table and adjustment knob.

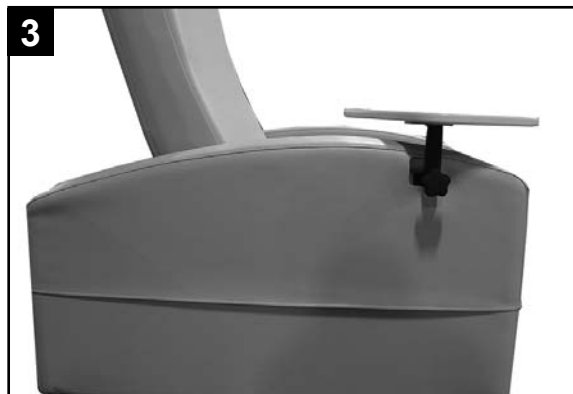

# **LIMITED WARRANTY**

## **TERMS**

Winco Mfg., LLC. warrants this product to be free of manufacturer's defects in material and workmanship, provided that the product is used according to normal operating conditions and proper maintenance intended by manufacturer. This warranty is available only to the original retail purchaser, is non-transferable and commences on the date of retail sale; proof of purchase required.

**LIMITED ONE YEAR COMPLETE PRODUCT WARRANTY:** Winco Mfg., LLC warrants the complete product for one (1) year.

At Winco Mfg., LLC. sole discretion, it may repair or replace any components freight free that are found to be defective during the first year. Winco Mfg., LLC. shall not be liable for any labor, or any other costs incurred as a result of or in conjunction with a warranty claim.

**LIMITED THREE YEAR WARRANTY:** Winco Mfg., LLC warrants recliner mechanisms, electrical components, vinyl panels and

Vinyl on upholstered parts [from cracking or delaminating] pressurized gas springs, and casters, when new, for a period of three (3) years. At Winco Mfg., LLC.'s Sole discretion, it may repair or replace components found to be defective. Winco Mfg., LLC. shall not be liable for any labor, shipping or any other costs incurred as a result of or in conjunction with a warranty claim.

**LIMITED LIFETIME CHAIR FRAME WARRANTY:** Winco Mfg., LLC warrants steel and aluminum base frames for all products,

when new, for the lifetime of the original purchaser. This limited warranty does not apply to paint/finish or any components attached to the frame such as; upholstery, foam, casters, mechanisms or related parts that are covered under above warranties. Winco Mfg., LLC shall not be liable for any labor, shipping or any other costs incurred as a result of or in conjunction with a warranty claim.

The purchaser's exclusive remedy under this warranty shall be limited to such repair or replacement of defective components at Winco Mfg., LLC sole discretion. For warranty service, contact the dealer that the product was originally purchased from or Winco Mfg., LLC directly.

## **EXCLUSIONS**

There are no other warranties, conditions, representations or guarantees, express or implied, made or intended by Winco Mfg., LLC and all other warranties, conditions, representations or guarantees including any warranties, conditions, representations or guarantees under any Sale of Good Act or Like legislation or statute is hereby expressly excluded. Any and all other implied warranties shall not extend beyond the duration of the express warranty. Liability for incidental or consequential damages is excluded to the extent permitted by law. Some states do not allow the exclusion or limitation so the above limitation or exclusion may not apply to you. This warranty gives you specific legal rights. In addition you may also have other rights, which vary from state to state.

## **GENERAL PROVISIONS**

No warranty herein contained or set out shall apply when damage or repair is caused by any of the following:

1. Damage in transit.
2. Accident, alteration, abuse or misuse of product i.e. exceeding weight capacities, applying inordinate pressure to footrest/leg rest, use in Mental Health facilities like institutions or, any unintended use of the product or use in unintended environments (i.e. outdoor, showers, MRI rooms, etc..)
3. Fire, water damage, theft, war, riot, hostility, acts of God.

***Examples contained in this list are not to be construed as all-inclusive.  
Design, appearance, parts, & construction are subject to change without notice.***



5516 Southwest 1st Lane  
Ocala, Florida 34474

1-800-237-3377/ 352-854-2929/ Fax: 352-854-9544

e-mail: [customerservice@wincomfg.com](mailto:customerservice@wincomfg.com)

[www.wincomfg.com](http://www.wincomfg.com)