

8570

Adjustable Headrest Treatment Table



CUSTOMER INSTRUCTIONS

PLEASE READ AND FAMILIARIZE YOURSELF WITH ALL INSTRUCTIONS BEFORE USING THIS PRODUCT. If you have trouble understanding these instructions contact your dealer or Winco customer support, (800) 237-3377 before attempting to use this product; otherwise injury may occur.

Winco assumes no responsibility for damage or injury caused by improper assembly, installation, use, or maintenance of these products.



1. READ AND FOLLOW ALL DIRECTIONS.
2. NEVER stand on table or INJURY MAY OCCUR.
3. NEVER allow user to sleep unattended on table or INJURY MAY OCCUR.
4. NEVER leave an unconscious user unattended on table or INJURY MAY OCCUR.
5. DO NOT use table for transporting user in or with ANY type of vehicle or trailer. Winco tables have not been tested or approved for use by an occupant in any type of vehicle or trailer.
6. NEVER leave an infant or child unattended on table or INJURY MAY OCCUR.
7. DO NOT put hands, feet or clothing into any openings when changing positions on table. Attendant MUST confirm that users hands and feet are safe while changing table positions or INJURY MAY RESULT.
8. Periodically, recheck tightness of all nuts, bolts and screws.
9. USE ONLY WINCO AUTHORIZED REPLACEMENT PARTS.
10. NEVER EXCEED the recommended **weight capacity of 400 lbs.** (181.4 kg).
11. SAVE THESE INSTRUCTIONS for future reference and training.

8570 Weight Capacity = 400 lbs. (181.4 kg.)

Tools needed:

- Phillips-head screwdriver
- 7/16" wrench (multiple size wrench included)
- Scissors or box cutter
- Cardboard shipping container (stays under table during assembly to protect upholstery)
- - If using powered tools always start hardware together by hand first & wear safety goggles

PLEASE READ AND FAMILIARIZE YOURSELF WITH ALL INSTRUCTIONS BEFORE PROCEEDING WITH ASSEMBLY

1. After checking your product for any shipping damage, table needs to be positioned upside-down on the cardboard shipping container, or protective surface. Remove all packing material and hardware that was secured for shipping. Cut tape on corners of box so that cardboard will lay flat while you assemble table.
2. You will find two **Welded Leg Assemblies**, one is attached to table frames and one is not. (*SEE parts guide*)
3. Pivot up the attached **Welded Leg Assembly**, (located on opposite end from the **Backrest**.)
4. Next, using a Phillips screwdriver, install (2) ¼-20 x 1-½" screws through the upper leg brace of the **Welded Leg Assembly** into (2) threaded holes in the wood bottom of **Upholstered Top**, as shown.
5. Install the unattached **Welded Leg Assembly**, on the end near the Backrest, as shown in parts guide. The cross brace needs to be positioned toward "backrest end" of table.
6. Install (2) ¼-20 x 2-½" screws and (2) ¼-20 nuts to fasten the unattached **Welded leg Assembly** to the **Table Frame Rails**. Heads of screws should face the sides of table, for better appearance. DO NOT FULLY TIGHTEN AT THIS TIME.

7. REPEAT Step 4
8. Place (4) ¼-20 x 1-½” screws in (4) holes through the lower leg brace of the **Welded Leg Assemblies** as shown. Heads of screws should face the ends of table, for better appearance.
9. Align **Welded Brace Assembly** (smoothest side up, when the table is in the upright position) with the (4) screws & holes on both **Welded Leg Assemblies** and fasten with (4) ¼-20 nuts.

If the **Optional Shelf** and **Reinforced Welded Brace Assembly** have been ordered; make sure that the smooth surface (side with tape) of the **Reinforced Welded Brace Assembly** will face upward (when the table is in the upright position) prior to fastening to **Welded Leg Assemblies**.

10. TIGHTEN ALL FASTENERS AT THIS TIME. (including the nuts & screws holding the Welded Leg Assemblies to the Table Frame Rails – see sketch).
11. Install the four (4) **Glide Feet** completely into the ends of the **Welded Leg Assemblies**.
12. The release handle for the adjustable back rest, can be installed on the left or the right table frame rail, depending on user preference.
 - Should you wish to change the release handle mounting location, simply remove the (2) Phillips head screws that attach the handle to the table frame rail and relocate the handle to the opposite table frame rail and re-install the screws.
13. Using at least two people, turn table upright.
14. Slide the **Paper Holder Rod** into the (2) **Paper Rod Clamps**.
15. Adjust one or two of the **Glide Feet** to level table.


Optional Shelf installation:

16. Double-sided tape has been installed on the Welded Brace Assembly.
17. Before removing paper on tape, “test fit” the shelf. Shelf should be approximately ¼” from edge of Lower Leg Brace.
18. Once you are certain of where you will position the shelf, peel off the paper to expose adhesive.
19. Carefully position shelf in place on adhesive.

*After peeling paper from adhesive tape:
Use damp sponge to moisten tape to allow shelf to be adjusted during install.
Once tape is dry shelf will stay in place.*

Operation of Table

- **To raise** the Upholstered Backrest, pull black lever on the Remote Release System.
- **To lower** the Backrest, pull the same black lever with one hand and apply pressure to the Backrest with your other hand until the desired position is achieved, then release the black lever.

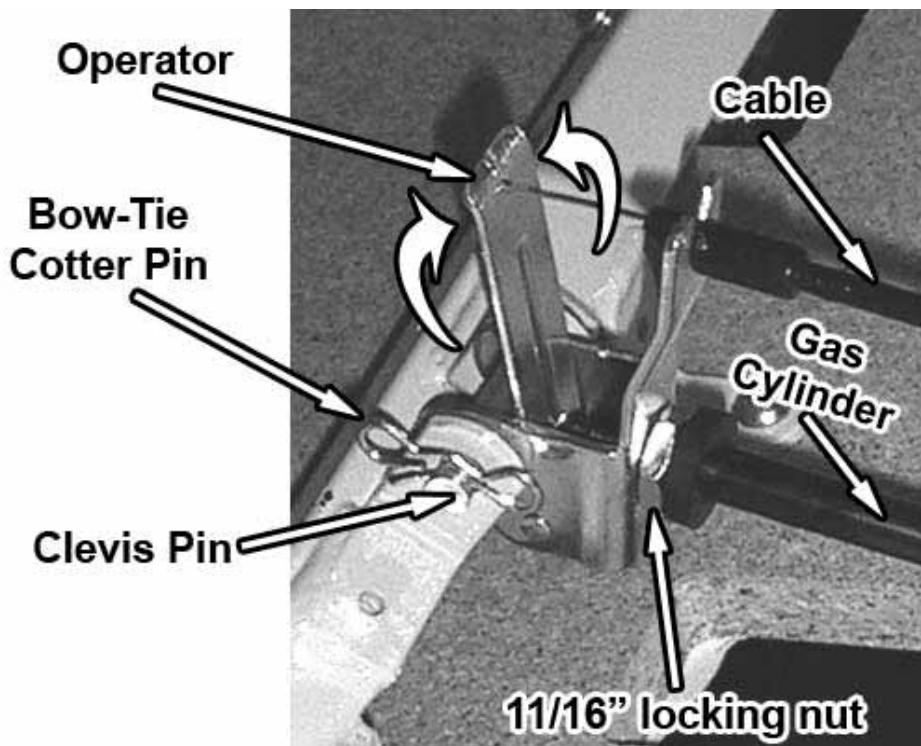

ADJUSTING GAS SPRING (cylinder):

1. Remove **Cable** from **Operator**.
2. Loosen **11/16"** nut on top of **Operator**.
3. Remove **Bow-Tie Cotter Pin** from **Clevis Pin**.
4. Remove **Clevis Pin**.
5. Turn the entire **Operator** one full revolution:
Clockwise if backrest does not move when handle is actuated OR Counter-clockwise if backrest will not stay in down position.

Finish the adjustment:

6. Re-install the **Clevis Pin**.
7. Re-install the **Bow-Tie Cotter Pin**.
8. Tighten the **11/16"** locking nut on top of **Operator**.
9. Re-install **cable**.
10. REPEAT these steps as needed until proper backrest action is achieved.



Vinyl Strap Adjustment – Should the vinyl strap on your Winco table need to be adjusted, unscrew the screw that currently holds either side of it in place. Pierce another hole in the strap to give desired tension and re-insert screw.

Wrinkles in Vinyl – Storage and shipping can cause wrinkles of the vinyl on your Winco table. This is normal. Please allow 1-2 weeks at room temperature (approx. 70 degrees) to allow for vinyl to expand back to original appearance.



Do not over tighten nuts and bolts. Periodically check all fasteners for tightness.

Parts Guide and List

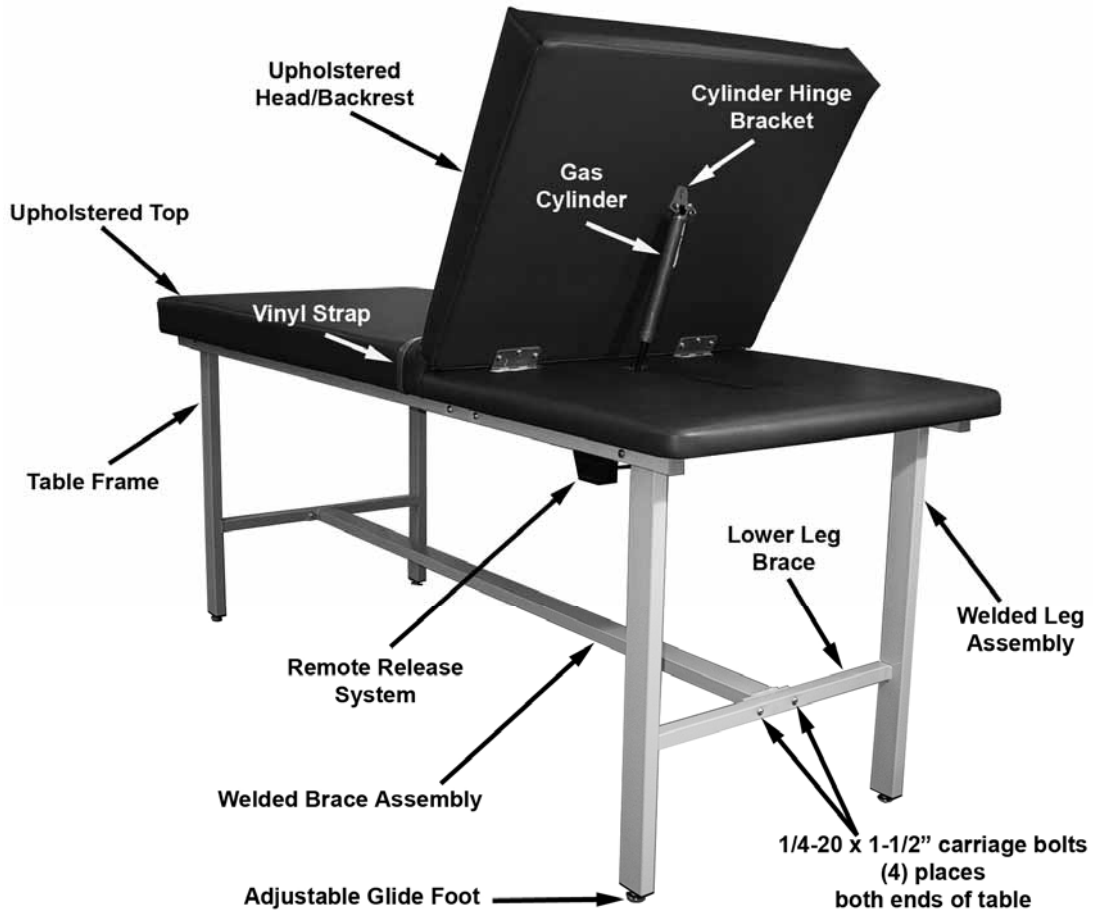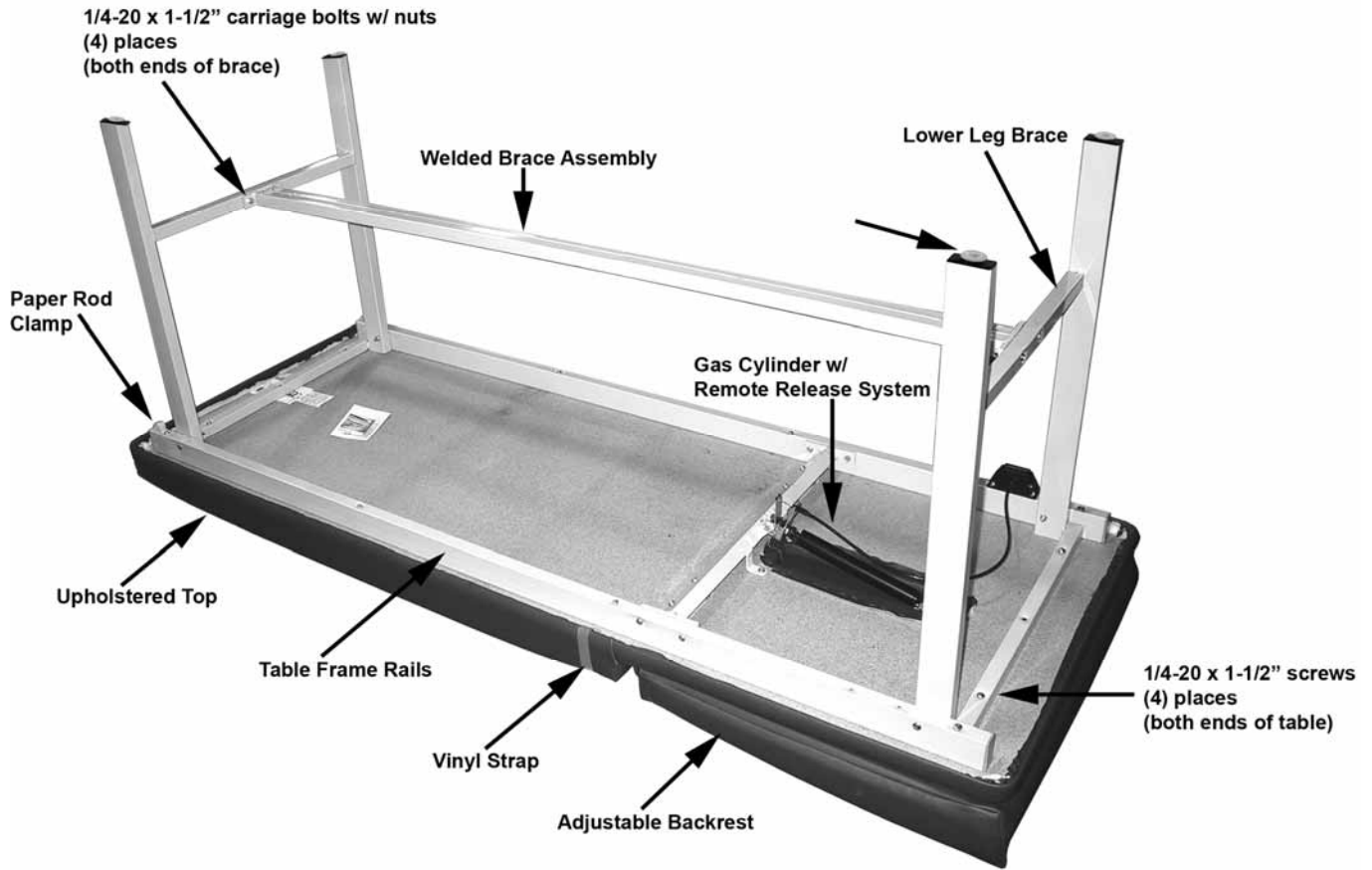

Parts Guide and List (cont'd)


<u>Items factory installed on table:</u>	<u>Qty.</u>	<u>Sizing chart – part number</u>
Table Frame Rails	2	222801
Paper Rod Clamp	2	501810
Cylinder Bracket Weldment	1	222830
Backrest Hinge	2	301250
Vinyl Strap (for paper)	1	271500
Gas Spring (cylinder) w/ Operator Assm. & Adj.	1	205225
Gas Cylinder	1	405225
Remote Release System	1	403450
Grommet, hinge (blk.)	2	700255
Grommet, slot 12 in. (blk.)	1	700260
Glide Insert Caps	4	701010
Upholstered Top	1	CALL FOR SPECIAL PART NUMBERS
Upholstered Head/Backrest	1	CALL FOR SPECIAL PART NUMBERS

Items included with table (not factory installed):

Welded Brace Assembly (used without shelf)	1	222720
Welded Brace Assembly (used with shelf)	1	222730
Shelf, plain	1	263960
Welded Leg Assembly (30" standard)	2	222710
Special Height Leg Assembly		CALL FOR PART NUMBERS
Glide Feet (leveling feet)	4	500350
Threaded insert caps (black)	4	500105
Paper Holder Rod	1	222705
Plastic Cap (white)	2	700650

¼ - 20 x 1½ carriage bolt
(4) leg assys. to cross brace

4

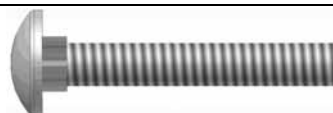

566401

¼-20 x 2½" Machine Screw
(1) leg assembly to table frame

2



566519

¼-20 x 1-½" Machine Screw
(4) leg assys. to table bottom

4



566506

¼-20 Nylon Lock Nut
(2) leg assy to table frame
(4) leg assys. to cross brace

6



566204

Other Quality Winco Medical Furnishings You May Be Interested In:



Clinical Recliners




Exam Tables



Convalescent (LTC) Recliners



Ring Binder Carts

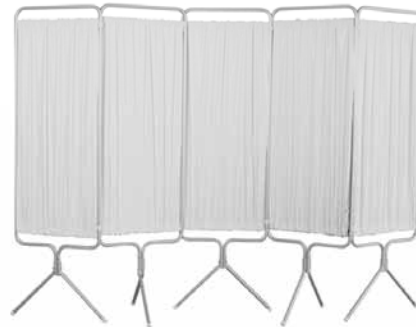

Privacy Screens



Phlebotomy Chairs



Specialty Chair Transfer Chairs



Stools



Room Chairs

Visit us online at www.wincomfg.com



ISO 9001-2008 CERTIFIED



LIMITED WARRANTY

• 5516 SW 1st Lane • Ocala, FL 34474
 • Phone (352) 854-2929 • (800) 237-3377
 • Fax (352) 854-9544
 • www.wincomfg.com
 • customerservice@wincomfg.com

TERMS

Winco Mfg., LLC. warrants this product to be free of manufacturer's defects in material and workmanship, provided that the product is used according to normal operating conditions and proper maintenance intended by manufacturer. This warranty is available only to the original retail purchaser, is non-transferable and commences on the date of retail sale; proof of purchase required.

Limited One Year Complete Product Warranty: Winco Mfg., LLC warrants the complete product for one (1) year. At Winco Mfg., LLC. sole discretion, it may repair or replace any components freight free that are found to be defective during the first year. *Winco Mfg., LLC. shall not be liable for any labor, or any other costs incurred as a result of or in conjunction with a warranty claim.*

Limited Three Year Warranty: Winco Mfg., LLC warrants the recliner mechanisms, electrical components, vinyl [from cracking or delaminating] pressurized gas springs, when new, for a period of three (3) years. At Winco Mfg., LLC. sole discretion, it may repair or replace components found to be defective. *Winco Mfg., LLC. shall not be liable for any labor, shipping or any other costs incurred as a result of or in conjunction with a warranty claim.*

Limited Lifetime Chair Frame Warranty: Winco Mfg., LLC warrants the steel chair frame, when new, for the lifetime of the original purchaser. This limited warranty does not apply to paint/finish or any components attached to the frame such as; upholstery, foam, casters, mechanisms or related parts that are covered under above warranties. *Winco Mfg., LLC shall not be liable for any labor, shipping or any other costs incurred as a result of or in conjunction with a warranty claim.*

The purchaser's exclusive remedy under this warranty shall be limited to such repair or replacement of defective components at Winco Mfg., LLC sole discretion. For warranty service, contact the dealer that the product was originally purchased from or Winco Mfg., LLC directly.

EXCLUSIONS

There are no other warranties, conditions, representations or guarantees, express or implied, made or intended by Winco Mfg., LLC and all other warranties, conditions, representations or guarantees including any warranties, conditions, representations or guarantees under any Sale of Good Act or Like legislation or statute is hereby expressly excluded. Any and all other implied warranties shall not extend beyond the duration of the express warranty. Liability for incidental or consequential damages is excluded to the extent permitted by law. Some states do not allow the exclusion or limitation so the above limitation or exclusion may not apply to you. This warranty gives you specific legal rights. In addition you may also have other rights, which vary from state to state.

GENERAL PROVISIONS

No warranty herein contained or set out shall apply when damage or repair is caused by any of the following:

1. Damage in transit.
2. Accident, alteration, abuse or misuse of product ie. exceeding weight capacities, applying inordinate pressure to footrest/leg rest, use in Mental Health facilities or like institutions.
3. Fire, water damage, theft, war, riot, hostility, acts of God.

Examples contained in this list are not to be construed as all-inclusive.

