



*The Winco Recliner Collection™*

# Augustine Treatment Recliner



We believe the Augustine Treatment Recliner can truly make a difference in patient care. The Augustine was developed to be 100% patient-centric; to be comfortable and help put the patient at ease. The result is an attractive and inviting, clinically preferred recliner, with unique solutions for enhanced clinical utility.

- Enhanced range of motion for optimum, patient selectable, positioning
- Spacious seat creates personal space and luxurious comfort
- Unmatched stability coupled with thoughtful features to ease chair movement and cleaning
- Styling and convenience enhanced with nearly invisible concealed casters

Developed in conjunction with Oncology clinics, Augustine is also ideally suited for use in patient rooms, recovery, and Long Term Care. Augustine demonstrates Winco's commitment to developing products that meet the needs of the people that live in and touch the recliners on a daily basis.

## SERENITY: PERSONAL SPACE AND INDEPENDENCE

Often the medical treatment environment unavoidably impinges upon the patient's personal space. Augustine addresses this with a welcoming and spacious patient seating area.

- Low seat height to ease entry and exit.
- Room to stretch out or curl-up and relax
- Easy to recline which allows desired positioning without assistance.

Patients were clear when asked their needs; they wanted a chair where they could feel at ease, have some independence and a relaxing experience. Augustine is designed to mimic furniture that you would be proud to have in your own home.

Augustine is designed to mimic furniture that you would be proud to have in your own home.



*The Central Front Braking option (CFB).*

## CLINICIAN INPUT: DESIGN ENHANCEMENT

The quick release seat is an Augustine clinician-driven solution. Facility managers requested an easier way to thoroughly clean the recliner. Once removed, full access to the inner sides of the recliner, backrest and seat allows for fast and thorough cleaning.

Simply pull-up at the front and slide forward. There is no need for any tools, and no risk of lost parts. The seat is securely tethered so it won't get separated from the recliner. To replace simply line up the tabs, slide into place and push down – done!

**REMOVABLE SEAT:** There is no need for any tools or risk of lost parts.



## CLINICAL UTILITY: PRACTICAL INNOVATION

Augustine has the clinical features required in a high use medical environment – wheels for positioning, tested durability with "breakaway" mechanism and ease of cleaning. These chairs bring dignity, comfort and grace back to the clinical setting while forfeiting none of the utility.

- Standard rear central caster lock offers impressive stability.
- Efficient thorough cleaning with removable seat, no tools required, and no lost parts.
- Maximize clinic space with exclusive new Central Front Brake Option (CFB).

CFB allows effortless, secure braking of the chair from the front, with no obtrusive levers or mechanisms to catch on or impede staff, patients or visitors.






The Winco Recliner Collection™

# Augustine Treatment Recliner

## ANATOMY OF THE AUGUSTINE TREATMENT RECLINER

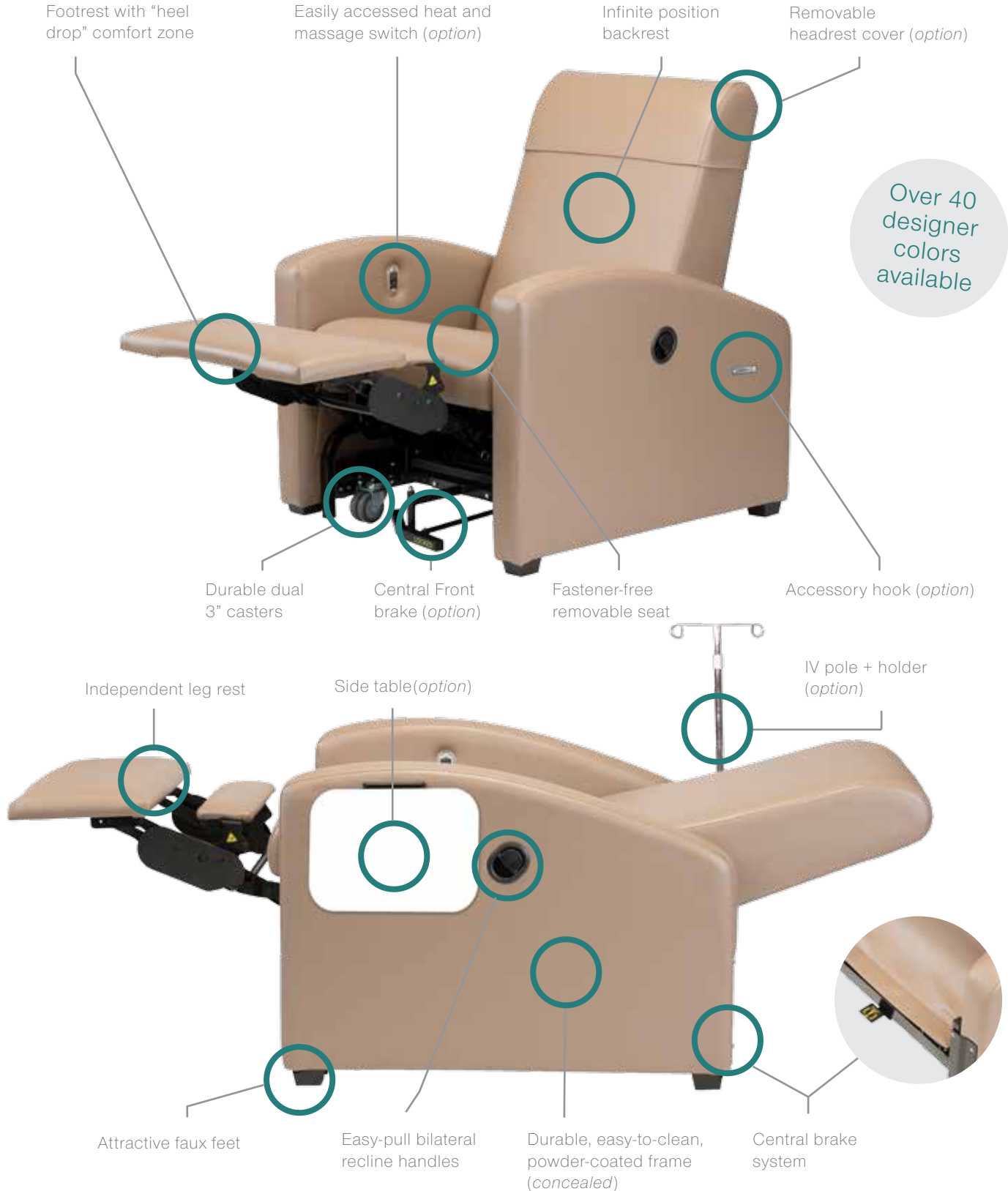

## SPECIFICATIONS

Seat Height (From Floor)	19.75" (50.2 cm)
Seat Width (Between Arms)	26" (66 cm)
Seat Depth	20.75" (52.7 cm)
Overall Height	45.75" (116.2 cm)
Overall Width*	34.5" (87.6 cm)
Back Height (Above Seat)	28.5" (72.4 cm)
Weight	181 lbs (82.1 kg)
Weight Capacity	500 lbs (226.8 kg)
Footprint (upright)	34.5" x 40" (87.6 cm x 101.6 cm)
Footprint (reclined)	34.5" x 76" (87.6 cm x 193 cm)
Seat-to-Back Angle Range	75° (degrees)

\* Add 2.5" to chair width for each side table in down position.  
Add 11.0" to chair width for each raised side table.

Over 40  
designer  
colors  
available



## FEATURES

MODEL	CASTER	ARM STYLE	STYLE	POSITIONING
5001	Front: 3" Fixed Rear: 3" Swivel Central Locking	Fixed	Standard	Infinite

## OPTIONS

ACCESSORIES			TABLES		COMFORT	BRAKING
IV Pole	Accessory Hook	Flame Rating	Side Table	Pivot Table	Heat Only or Heat & Massage	Central Front Brake
(L, R or Both)	(L, R or Both)	TB133	(L, R or Both)	(L, R or Both)	•	•

## RECOMMENDED APPLICATIONS

ONCOLOGY & INFUSION	ORTHOPEDICS	PATIENT ROOM	SAME-DAY SURGERY	LONG-TERM CARE
•	○	•	○	○

• Optimal Application ○ Suitable Application

**WARRANTY** Contact Winco or visit our website for warranty details.

**VINYL** Contact Winco or visit our website for available vinyl selections.

**NOTE** Winco reserves the right to change specifications without notice.

**CERTIFICATIONS** **1)** The Augustine Treatment Recliner meets ANSI/BIFMA X5.4 for reliability testing **2)** Winco foam meets California Technical Bulletin 117, Flammability Test Procedure for Resilient Filling Materials used in Upholstered Furniture **3)** TB133 Option: Meets flammability standards of California Technical Bulletin 133, Flammability Test Procedure for Seating Furniture for Use in Public Occupancies and ASTM E1537, Standard Test Method for Fire Testing of Upholstered Furniture **4)** Heat and Massage, Heat Only Options: Meets electrical safety standards, UL 1647 and CAN/CSA-C22.2 No.68-92

\* Patent Pending



## Contact Us

To arrange a demonstration or to speak with one of our associates, please feel free to contact us at:  
800-237-3377 | WincoMfg.com | 5516 SW 1st Lane, Ocala, FL 34474 | Info@WincoMfg.com